


Immersion tube

This document describes how to mount the immersion tube on a HAAKE™ Viscotester™ iQ rheometer.

Mounting

The immersion tube is designed for measuring a sample in a larg(er) (original) container using the well defined CC25 DIN measuring geometry.

Figure 1. Immersion tube

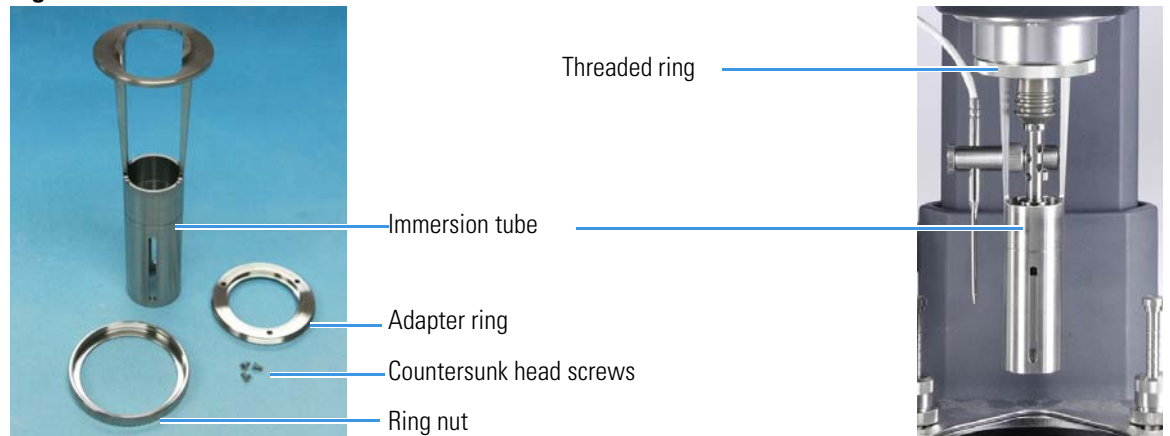

❖ To mount the immersion tube

1. Mount the adapter ring to the flange at the bottom side of the instrument head using the 3 countersunk screws.
2. Attach the immersion tube to the adapter ring.
3. Fasten the immersion tube to the instrument head by screwing the ring nut to the adapter ring.

❖ To dismount the immersion tube

1. Loosen the ring nut from the adapter ring.
2. Remove the immersion tube from the instrument head.
3. The adapter ring must not be removed from the instrument.